



DOC023.53.90050

TSS Portable handheld measurement instrument for turbidity/solids

User Manual

06/2012, Edition 3

Table of Contents

Section 1 Specifications	5
Section 2 General Information	7
2.1 Safety information	7
2.1.1 Hazard notices in this manual	7
2.1.2 Precautionary labels	7
2.2 Overview of product	8
2.3 Measurement instrument	8
2.4 Measuring principle	8
2.5 Probe	8
Section 3 Installation	9
3.1 Product contents	9
3.2 Powered via rechargeable batteries	10
3.2.1 Inserting the rechargeable batteries	10
3.2.2 Install the charger	11
3.2.3 Charging the rechargeable batteries	12
3.3 Connecting the probe	13
3.4 Switching the instrument on and off	14
3.4.1 Register a new probe	14
Section 4 Start Up	15
4.1 System start overview	15
4.2 User interface and navigation	15
4.2.1 Keypad	15
4.2.2 Display	16
4.2.2.1 Querying internal information	17
4.3 Select language	17
4.4 Setting the date and time	17
4.5 Setting the display	18
4.6 Setting the units	18
4.6.1 Sludge blanket level	18
4.7 Calibration	19
4.8 Defining the measurement points	19
4.8.1 Overview of calibration curves	19
4.8.2 Measurements	19
4.8.2.1 Practical example for solids measurement	19
4.8.2.2 Practical example for turbidity measurement	21
4.9 Setting the integration time	21

Section 5 Operation	23
5.1 Calibration	23
5.1.1 Important notes for calibration	23
5.1.2 Calibrating	24
5.1.3 Manual calibration value correction	26
5.1.4 Delete calibration point	26
5.1.5 Reset to default calibration	26
5.2 Measure	26
5.2.1 Selecting the calibration curve	27
5.2.2 Start the measurement	27
5.2.3 Interrupt and resume the measurement	27
5.2.4 Stop measurement	27
5.3 Show the data in the display	27
5.4 Delete data for a calibration curve	27
5.5 Delete stored data for all calibration curves	28
Section 6 Maintenance	29
Section 7 Troubleshooting	31
7.1 Error messages	31
7.2 Informative messages	31
Section 8 Replacement parts and accessories	33
8.1 Replacement Parts	33
Section 9 Contact information	35
Section 10 Limited warranty	37

Section 1 Specifications

These are subject to change without notice.

Performance specifications		
Wavelength	860 nm	
Parameter	Turbidity	Solids (dry matter)
Measurement method	Combined multiple-beam alternating light technique with IR diode system and beam focus	
	2-channel 90° scattered light measurement in accordance with DIN EN 27027 / ISO 7027; additional measurement value verification through six-channel multiple-angle measurement	Modified absorption measurement Six-channel multiple-angle measurement
Measuring range	0.001–9999 FNU (NTU)	0.001–400 g/L (upper limit depends on matrix)
Resolution	0.001 at 0–0.999 FNU 0.01 at 1–9.99 FNU 0.1 at 10–99.9 FNU 1 at >100 FNU	0.001 at 0–0.999 g/L 0.01 at 1–9.99 g/L 0.1 at 10–99.9 g/L 1 at >100 g/L
Precision	Measuring range: 0.001–9999 FNU Accuracy of the measurement value: < 3% or +/- 0.02 FNU (whichever is greater)	Measuring range: 0.001–400 g/L Accuracy of the measurement value: < 4% or +/- 0.001 g/L (whichever is greater)
Reproducibility	<4% of measurement value	<5% of measurement value
Units	FNU, NTU, EBC	ppm, mg/L, g/L, %
Calibration	1 calibration curve (factory calibration)	4 calibration curves
Sample temperature	0–60 °C (32–140 °F), up to 80 °C (176 °F) briefly	
Pressure range	Max. 10 bar	
Display	LCD, alphanumeric, 4 rows of 16 characters	
Input	6 membrane keys, menu with quick access to important functions	
Power supply (rechargeable batteries)	6 rechargeable NiMH batteries (recommended: 1.2 V/min. 1800 mAH)	
Power consumption	Approx. 60 mA	
Data log	Up to 290 measurement values	
Interface	RS 485	
Probe material	Stainless steel, sapphire	
Cable	10 m (33 ft), PUR, Ø 8.3 mm (0.33 in.); S-2000 connector, 6-pin	
Protection class	Probe: IP68 Control unit: IP 55	
Size	Probe: Ø 40 mm (1.57 in.), length = 29 cm (11.42 in.) Control unit: 11 x 23 x 4 cm (4.33 x 9 x 1.57 in.)	
Weight	Probe: 1600 g (3.53 lb) Control unit: 560 g (1.23 lb)	
Warranty	2 years	

Section 2 General Information

2.1 Safety information

Please read this entire manual carefully before unpacking, setting up, or operating this equipment. Pay attention to all danger and caution notices. Failure to do so could result in serious injury to the operator or damage to the equipment.

To prevent damage to or impairment of the device's protection equipment, the device may only be used or installed as described in this manual.

2.1.1 Hazard notices in this manual

DANGER

Indicates a potentially or imminently hazardous situation that, if not avoided, can result in death or serious injury.

WARNING

Indicates a potentially or imminently dangerous situation that, if it is not avoided, can lead to death or to serious injuries.

CAUTION

Indicates a potentially hazardous situation that may result in minor or moderate injury if it is not avoided.




NOTICE

Indicates a situation that, if it is not avoided, can lead to damage to the device. Information that requires special emphasis.

Note: Information that supplements points in the main text.

2.1.2 Precautionary labels

Read all labels and tags attached to the instrument. Personal injury or damage to the instrument could occur if not observed.

	This symbol, if noted on the instrument, references the instruction manual for operation and/or safety information.
	This symbol may be found on an enclosure or barrier within the product and indicates a risk of electric shock and/or death by electrocution.
	Electrical equipment marked with this symbol may not be disposed of in European domestic or public disposal systems after 12 August 2005. In conformity with European local and national regulations (EU Directive 2002/96/EC), European electrical equipment users must now return old or end-of life equipment to the manufacturer for disposal at no charge to the user. Note: Please contact the equipment manufacturer or supplier for instructions on how to return end-of-life equipment, manufacturer-supplied electrical accessories, and all auxiliary items for recycling or proper disposal.

2.2 Overview of product

The TSS Portable is a handheld measurement instrument for the analytical determination of turbidity and solids in aqueous media.

2.3 Measurement instrument

The instrument stores the recorded data under the corresponding calibration curve. Four calibration curves for solids (C-DS1, C-DS2, C-DS3, C-DS4) and one calibration curve for turbidity (C-TU) are available for selection.

For solids measurements, a specific calibration must be assigned to each measurement point ([section 5.1, page 23](#)).

All measurement values are saved with details of the selected calibration curve, the measurement value, the homogeneity, the date and the time.

Various individual parameters for input, signal processing and output can be set in the menu ([Section 4, page 15](#)).

2.4 Measuring principle

The measuring principle is based on a combined infrared absorption stray light process, which determines the lowest turbidity value according to DIN EN 27027 just as precisely and continuously as the high sludge content. In so doing, the light scattered sideways by the turbidity particles is measured at an angle of 90°. In the case of solid material, the measurement occurs at an angle of 90° and 120°.

2.5 Probe

The probe contains sensitive optical and electronic components. Care must therefore be taken to ensure that it is not subjected to any hard mechanical impacts. The inside of the probe, and of the display unit, does not contain any components that can be serviced by the user.

Section 3 Installation

⚠ WARNING

Risk of falling. Measurements are taken directly on site with this instrument. In order to avoid falling into the basin, all local safety regulations with respect to roping up and wearing appropriate protective clothing and shoes must be observed.

Risk of injury. Do not secure the instrument to your body. The measuring probe may catch unintentionally on a slide or stirrer and the connected probe cable may put the user at risk.

3.1 Product contents

After taking the instrument out of the case, check each part for damage. All the components on the enclosed list must be present. If any parts are missing or damaged, please contact the manufacturer or the distributor.

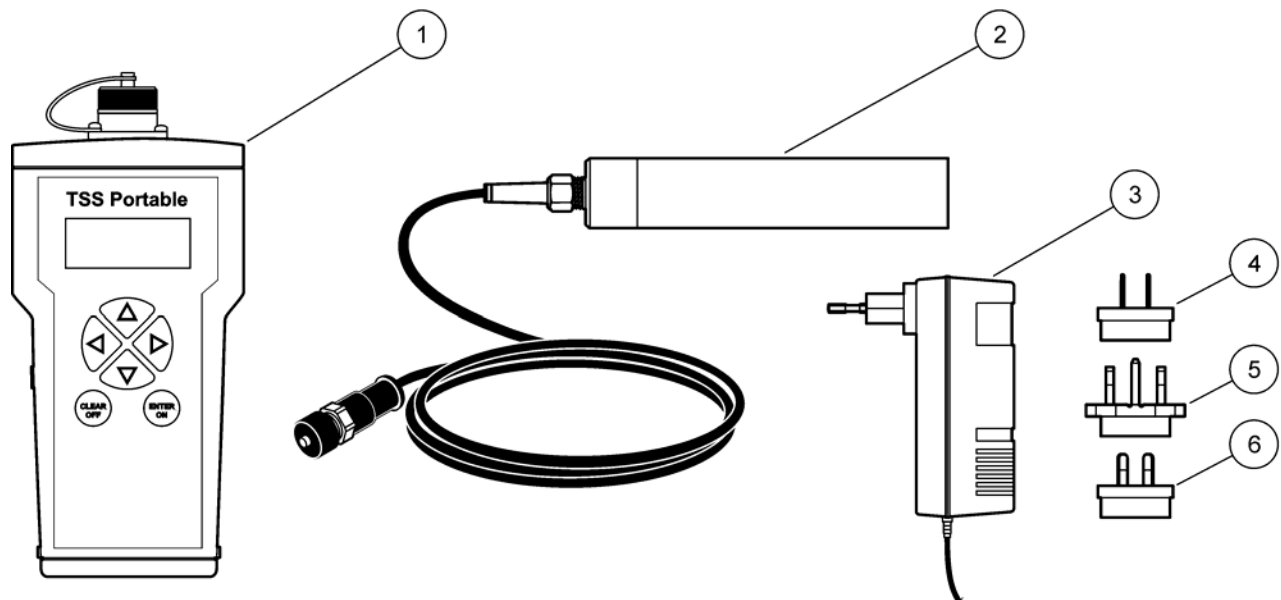


Figure 1 Product contents

1	TSS Portable handheld instrument	4	USA adapter plug
2	Probe with 10 m (32.8 ft) cable	5	UK adapter plug
3	Charger with EU adapter plug installed	6	AUS/China adapter plug

3.2 Powered via rechargeable batteries

⚠ WARNING

Risk of fire and explosion. Only use rechargeable NiMH batteries and ensure that the rechargeable NiMH batteries are inserted correctly in the battery compartment. Incorrect insertion of the rechargeable NiMH batteries may damage the instrument or lead to fire or explosions.

⚠ WARNING

Risk of fire and explosion. When using not rechargeable AA batteries, the not rechargeable batteries could explode if the instrument is connected to the charger. Only use rechargeable NiMH batteries from the manufacturer of the instrument.

The instrument is operated by six rechargeable NiMH batteries ([Figure 2](#)).

NOTICE

Remove the rechargeable NiMH batteries if the instrument is not to be used for a period of time. Only use rechargeable NiMH batteries from the manufacturer of the instrument.

3.2.1 Inserting the rechargeable batteries

1. Press in both catches at the same time and release the cover.
2. Take the battery holder out of the battery compartment ([Figure 2](#)).
3. Insert the delivered six rechargeable NiMH batteries into the battery holder. Take note of the markings for the polarity.

***Note:** Make sure that the rechargeable batteries have been inserted into the battery compartment correctly.*

4. Connect the battery holder to the battery clip.
5. Push the battery holder into the battery compartment.
6. Push the cover onto the battery compartment until it audibly engages in the catches.

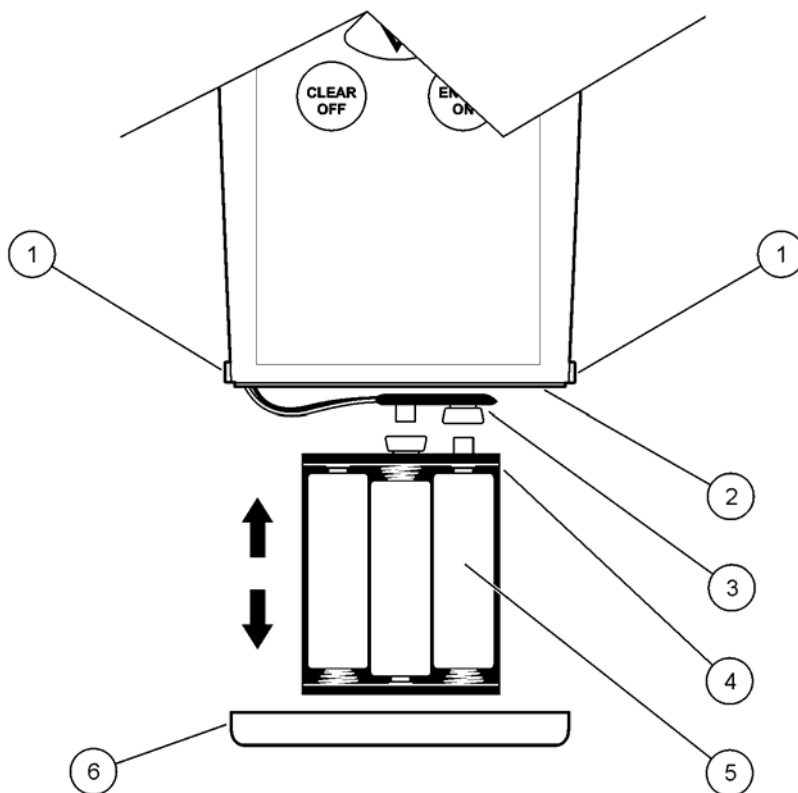


Figure 2 Inserting the rechargeable batteries

1	Catch	4	Battery holder
2	Battery compartment	5	rechargeable Batteries
3	Battery clip	6	Cap

3.2.2 Install the charger

⚠ WARNING

Risk of fire and explosion.
 Only use the battery charger LZY607 to charge the rechargeable NiMH batteries.

1. Push the slide switch (item 2, Figure 3) on the back of the charger up and remove the adapter (item 3, Figure 3).
2. Fix the required adapter to the charger, until it audibly locks into place.
3. Insert the charging plug (item 3, Figure 4) into the charging socket (item 4, Figure 4) of the instrument.
4. Connect the power (100–240 V~/50–60 Hz).

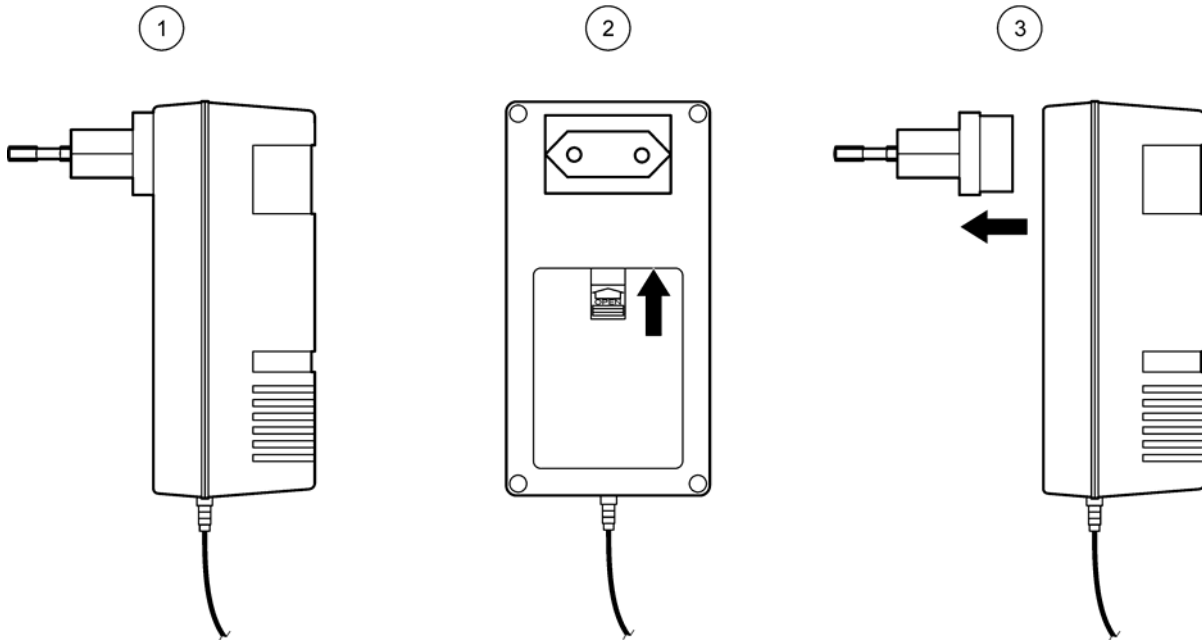


Figure 3 Replacing the adapter

1	Charger	3	Removing the adapter
2	Slide switch		

3.2.3 Charging the rechargeable batteries

Note: The instrument must be switched off while charging so that the rechargeable batteries can be charged.

Note: When using the instrument for the first time, charge for at least three hours. The rechargeable batteries are fully charged when the LED on the charger is green.

1. Connect the charger to the power supply source and then connect to the instrument (Figure 4). The rechargeable batteries are charged.
2. Once the rechargeable batteries are fully charged, disconnect the instrument from the charger.

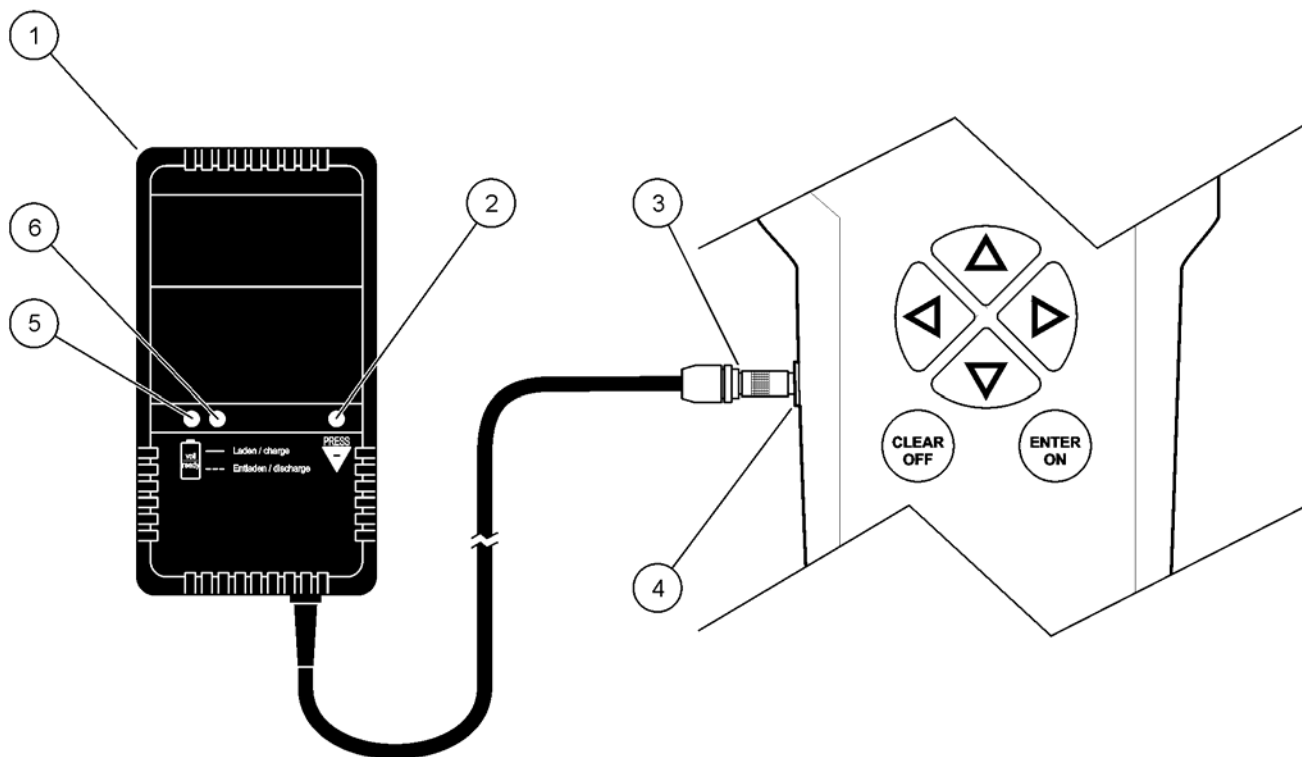


Figure 4 Charging the rechargeable batteries

1	Charger	4	Charging socket
2	Discharge button (yellow)	5	LED lamp (green => rechargeable batteries fully charged)
3	Charging plug	6	LED lamp (red => rechargeable batteries currently charging)

3.3 Connecting the probe

Remove the protective cap and connect the probe plug to the instrument ([Figure 5](#)).

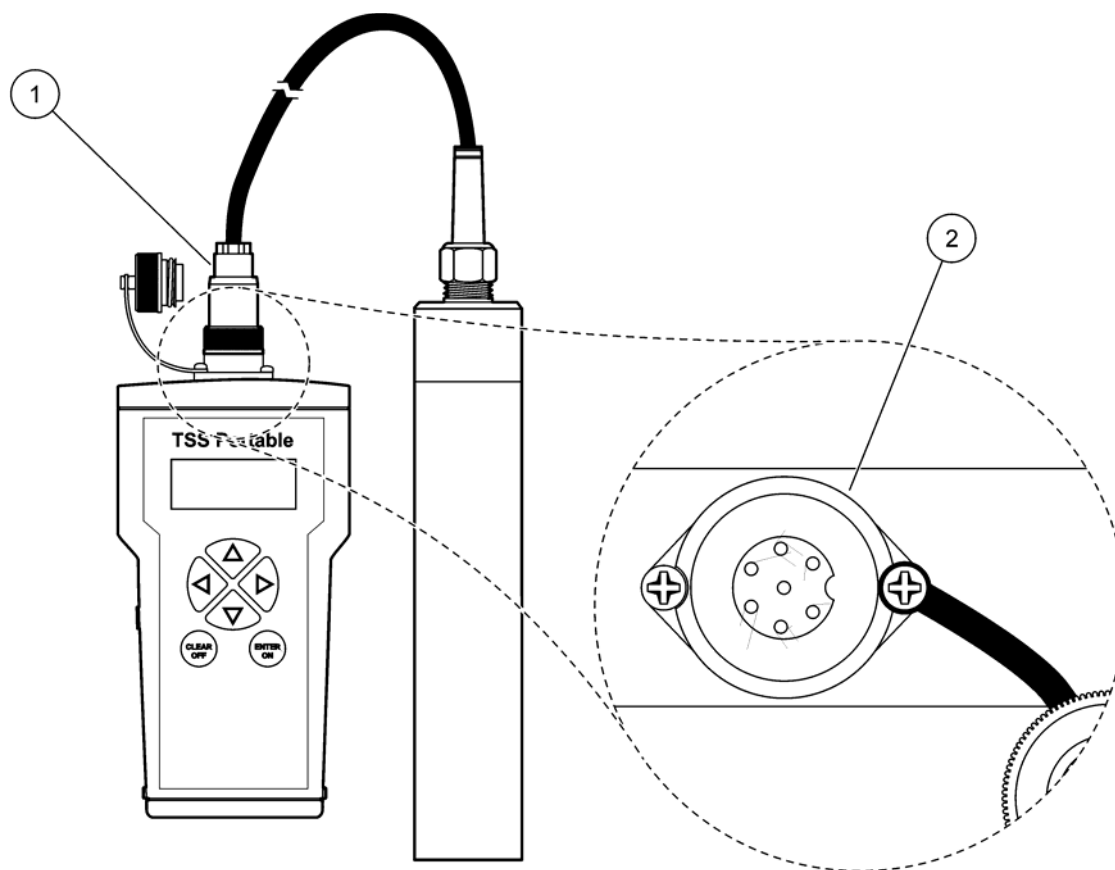


Figure 5 Connecting the probe plug to the instrument

1 Probe plug	2 Probe connection
--------------	--------------------

3.4 Switching the instrument on and off

1. Fully charge the rechargeable batteries (see [section 3.2.3](#)).
Note: The rechargeable batteries are fully charged when the LED on the charger is green.
2. Press **ENTER/ON** for two seconds to turn on the instrument.
Note: If the instrument does not turn on, check the position of the rechargeable batteries.
3. Press **CLEAR/OFF** for two seconds to turn the instrument off again.

3.4.1 Register a new probe

If a new probe is connected during initial installation, the error message "Incorrect probe" is displayed.

1. Select **MENU** and confirm with **ENTER/ON**.
2. Select **SYSTEM** and confirm with **ENTER/ON**.
3. Enter code **390** using the navigation keys (see [section 4.2, page 15](#)) and confirm with **ENTER/ON**.
4. Select **CLEAR/OFF** twice to return to the main menu.

Note: The probe is now registered and the TS curves can be calibrated (see [section 5.1, page 23](#)).

Section 4 Start Up

4.1 System start overview

1. Switch on the measurement instrument ([section 3.4, page 14](#)).
2. Connect the probe to the measurement instrument ([section 3.3, page 13](#)).
3. Edit the time/date, display, units, language and integration time ([section 4.3, page 17](#)).
4. Calibration of the solids calibration curve in line with calibration instructions ([section 4.7, page 19](#)).
5. Select desired calibration curve ([section 4.8, page 19](#)).
6. Show the data in the display ([section 5.3, page 27](#)).

4.2 User interface and navigation

4.2.1 Keypad

Figure 6 shows the navigation keys on the instrument.

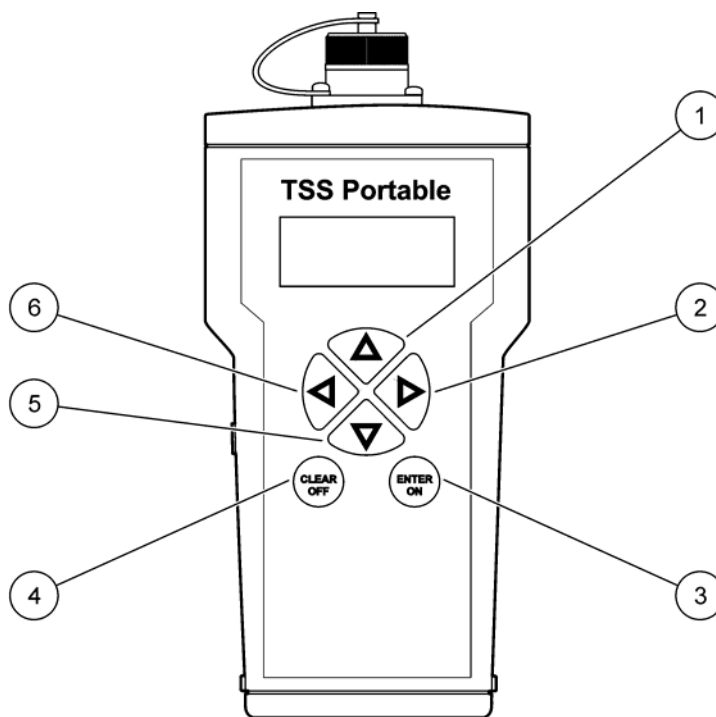


Figure 6 Description of the keypad

1	Navigation key UP	4	CLEAR/OFF
2	Navigation key RIGHT	5	Navigation key DOWN
3	ENTER/ON	6	Navigation key LEFT

4.2.2 Display

In its basic setting, the display shows the main measurement variables and the date/time (Figure 7). It can be set in a user-defined manner, see section 4.5, page 18.

Table 1 Navigating in the menu

Navigation key	Description
Navigation key RIGHT/LEFT	Navigation occurs with the RIGHT/LEFT navigation keys.
Navigation key UP/DOWN	Navigation occurs with the UP/DOWN navigation keys.
ENTER/ON	<ul style="list-style-type: none"> – Confirm selection and entry with ENTER/ON. – Switch the instrument on with ENTER/ON.
CLEAR/OFF	<ul style="list-style-type: none"> – Abort actions with CLEAR/OFF. – Press CLEAR/OFF to return to the previous menu option or abort entries. – Switch the instrument off with CLEAR/OFF.
Vertical arrows	The vertical arrows in the display show that there are additional menu options or measurement values above or below.
Cursor	<ul style="list-style-type: none"> – The cursor points to the selected main menu option. Navigation occurs with the RIGHT/LEFT navigation keys. – The active decimal place is displayed.
> (arrow)	The current selection is displayed

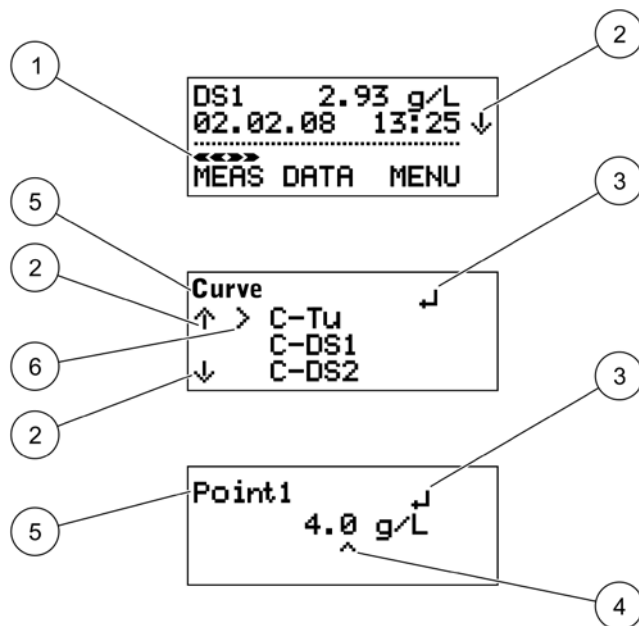


Figure 7 Display of results from the instrument with connected probe

1	Cursor	4	Current decimal place
2	Vertical arrows	5	Menu option
3	Input mode	6	Active selection

4.2.2.1 Querying internal information

Information about the probe and instrument can be queried directly from the main menu using the navigation keys **UP/DOWN**.

The following information is displayed:

- a. Measurement value (homogeneity in %)
- b. Battery charge level
- c. Diagnosis parameters, such as:
 - Probe serial number (instrument information)
 - Probe software version (instrument information)
 - Probe number Start XX (probe diagnosis data)
 - DATA rem. cap. XXX (Remaining capacity of the measurement points)
 - Fault no. XX
 - Serial no. XXXXX (instrument designation)
 - Software version XXXXX (instrument designation)

4.3 Select language

The instrument supports the languages German and English. The instrument functions in the selected language until the option is changed.

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **SPRACHE/LANGUAGE** and confirm with **ENTER/ON**.
3. Select **DEUT./GERMAN** for German or **ENGL./ENGLISH** for English and confirm with **ENTER/ON**.
4. Press **CLEAR/OFF** to return to the main menu.

4.4 Setting the date and time

Set the time and date display:

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **TIME/DATE**, confirm with **ENTER/ON**.
3. Select **YEAR**, **MONTH**, **DAY**, **HOUR**, **MINUTE** or **SECOND**, confirm with **ENTER/ON**.
4. Change the displayed values with the navigation keys.
5. Confirm each setting with **ENTER/ON**.
6. Press **CLEAR/OFF** to return to the **MENU**.

*Note: Press **CLEAR/OFF** again to return to the main menu.*

4.5 Setting the display

It is possible to custom-define the first two rows of the display. In the standard setting, the first line shows the main measured variable and the second line shows the date/time.

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **DISPLAY**, confirm with **ENTER/ON**.
3. Select **ROW 1**, confirm with **ENTER/ON**.
4. Select option for **ROW 1**, confirm with **ENTER/ON**.
 - Measurement value
 - Homogeneous
 - Battery
 - Time
5. Select **ROW 2**, confirm with **ENTER/ON**.
6. Select option for **ROW 2**, confirm with **ENTER/ON**.
7. Press **CLEAR/OFF** to return to the **MENU**.

*Note: Press **CLEAR/OFF** again to return to the main menu.*

4.6 Setting the units

The following measurement units are available for turbidity and solids:

Turbidity units (TU): FNU, NTU, EBC

Solids units (DS1 to DS4): ppm, mg/L, g/L, %

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **UNITS**, confirm with **ENTER/ON**.
3. Select the desired calibration curve (C-TU, C-DS1, C-DS2, C-DS3 or C-DS4) and confirm with **ENTER/ON**.
4. Select the desired measurement unit and confirm with **ENTER/ON**.
5. Press **CLEAR/OFF** to return to the **MENU**.

*Note: Press **CLEAR/OFF** again to return to the main menu.*

4.6.1 Sludge blanket level

To determine the sludge blanket level, set the units to mg/L, g/L or %. Select an integration time of 0 seconds (refer to [section 4.9, page 21](#)) to increase the response time. Calibrate the measurement system and slowly lower the sensor into the clarifier. Once the sludge blanket level is reached, the solids concentration will increase significantly. The distance of the sludge blanket from the surface of the water can then be read off from the meter markings on the probe cable.

4.7 Calibration

For turbidity measurement, a standard curve C-TU is already stored in the probe. Calibration is not required.

However, for solids measurement, calibration is required in order to set the rough measurement signals to a calibrated display (refer to [section 5.1, page 23](#)). A precise solids measurement is not possible without calibration. There are four calibration curves available: C-DS1, C-DS2, C-DS3, C-DS4. These curves can be allocated to the individual measurement points.

4.8 Defining the measurement points

4.8.1 Overview of calibration curves

- C-TU = Turbidity curve (standard curve)
- C-DS1 = solids curve 1 (calibration curve)
- C-DS2 = solids curve 2 (calibration curve)
- C-DS3 = solids curve 3 (calibration curve)
- C-DS4 = solids curve 4 (calibration curve)

4.8.2 Measurements

4.8.2.1 Practical example for solids measurement

Note: The type and composition of the solid particles of a medium can differ greatly. For this reason, it is not possible to define a general standard, and different measurement media must be calibrated specifically in each case.

For the solids measurement, calibration is required. There are four curves available: C-DS1, C-DS2, C-DS3, C-DS4. Each curve can be assigned to any measurement point. In the example, curve C-DS1 is selected.

1. Connect the probe to the instrument.
2. To switch on the instrument, press **ENTER/ON** for 2 seconds.

Select calibration curve

3. Select **MENU**, confirm with **ENTER/ON**.
4. Select **CALIBRATE**, confirm with **ENTER/ON**.
5. Select **CURVE C-TU** and confirm with **ENTER/ON**.

Note: Curve C-TU is the default setting.

6. Select calibration curve **C-DS1** with the navigation keys **UP/DOWN** and confirm with **ENTER/ON**.

The instrument reads the probe data and info text appears on the display.

Calibrate curve C-DS1:

7. Select **MEMORY** and confirm with **ENTER/ON**.
8. Lower the probe into a container that contains a homogeneous sample.
9. Select **POINT 1** and confirm with **ENTER/ON** while stirring the measuring medium with the probe.

The distance between the probe head and the walls and base of the container must always be more than 70 mm (2.76 in.).

It takes 5 to 20 seconds to record the calibration point; then the selection menu for calibration appears. While recording, the display shows the note "Memory".

10. Take the probe out of the container and clean it.
11. Determine the solids content of the sample in the laboratory.
12. Select ***POINT 1** and confirm with **ENTER/ON**.
13. Enter the laboratory value with the navigation keys and confirm with **ENTER/ON**.
 - Navigation key **LEFT/RIGHT**: Jump to next/previous decimal place
 - Navigation key **UP/DOWN**: Change number

The laboratory value is saved. The main menu is shown.

Measure:

14. Select **MEASURE**, confirm with **ENTER/ON**.
15. Select **START** and confirm with **ENTER/ON**.

Note: Data is automatically saved every minute, max. 290 measurement points.

Stop measurement:

16. Select **MEASURE**, confirm with **ENTER/ON**.
17. Select **STOP**, confirm with **ENTER/ON**.

Delete the saved measurement points:

18. Select **DATA**, confirm with **ENTER/ON**.
19. Select **DELETE MEP** and confirm with **ENTER/ON**.
20. Select **YES** and confirm with **ENTER/ON**.

4.8.2.2 Practical example for turbidity measurement

A standard curve C-TU is stored for the turbidity measurement.

1. Connect probe to instrument.
2. To switch the instrument on, press **ENTER/ON** for 2 seconds.
3. Select **READ** and confirm with **ENTER/ON**.
4. Place probe in measuring medium.
5. Select **START** and confirm with **ENTER/ON**.

Note: The data is automatically saved every minute, max. 290 measurement points.

6. To stop the measurement, select **READ** and confirm with **ENTER/ON**.
7. Select **STOP**, confirm with **ENTER/ON**.

Delete the saved measurement points:

8. Select **DATA** and confirm with **ENTER/ON**.
9. Select **DELETE MEP** and confirm with **ENTER/ON**.
10. Select **YES** and confirm with **ENTER/ON**.

4.9 Setting the integration time

Note: The measurement values are averaged out over the integration time, which results in a smooth measurement signal. For the sludge blanket level measurement (refer to [section 4.6.1, page 18](#)), select the integration time 0 seconds.

Upon delivery, the default value is set to 30 seconds. The integration time can be set to between 0 and 1000 seconds.

1. Select **MENU** and confirm with **ENTER/ON**.
2. Select **EXTRAS** and confirm with **ENTER/ON**.
3. Select **INTEGRATION** and confirm with **ENTER/ON**.
4. Change the displayed value with the navigation keys.
5. Confirm each entry with **ENTER/ON**.
6. Press **CLEAR/OFF** to return to the **MENU**.

*Note: Press **CLEAR/OFF** again to return to the main menu.*

5.1 Calibration

Turbidity measurements do not need to be calibrated, as a standardized calibration curve according to ISO 7027 is stored in the instrument. It is possible to create a customer-specific calibration of the turbidity curve. If the standard curve has been altered, a star (*) appears before the measurement point.

Solids measurements must be calibrated on site. It is not possible to define a general standard for this as the type and composition of the solid particles can differ greatly.

5.1.1 Important notes for calibration

Calibration containers

In order to avoid reflections of the measurement beam on the base and walls of the container, a black—or at the very least dark-colored—container should be used for calibration.

The probe can also remain in the tank so long as a representative laboratory sample can be taken at this point.

Measured variable (homogeneity)

100% represents a homogeneous measurement medium. After approximately 10 seconds, this measured variable provides a representative statement about the homogeneity of the measurement medium. The homogeneity is displayed as a percentage value after an individual measurement (Figure 8).

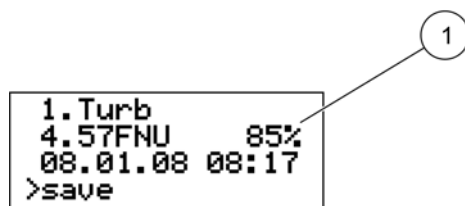


Figure 8 Homogeneity in %

1 Homogeneity

Replacing the probe

As the calibration is carried out for the measuring system (meter and probe), only the probe supplied may be used for the measurements. If a different probe is used, the equipment must be recalibrated.

Note: If you have changed the probe, the error message "Incorrect probe" is displayed. Register the new probe under **MENU > SYSTEM** with the code **390** (see [section 3.4.1, page 14](#)). This releases all the calibration curves for the solid matter measurement. The instrument then needs to be recalibrated.

5.1.2 Calibrating

The instrument can save one turbidity curve and up to four solids calibration curves. This makes measurements in media with different qualities possible. Each measurement point can have one of the saved calibration curves assigned to it individually.

1-point calibration

In order to define a calibration curve, it is generally sufficient to use a single calibration point in the probe's measuring range. This should lie in the upper third of the expected measuring range.

Note: *If the measurement value is below or above the range defined by a point, the instrument displays an error message: Calibration insufficient +/- . This means an additional calibration point is required.*

Multiple-point calibration

For broad measuring ranges, a second calibration point must be recorded:

- First calibration point (lower value) at as low a concentration as possible.
- Second calibration point (upper value) at as high a concentration as possible.

Up to three calibration points per curve can be recorded. To record several calibration points, the sample is diluted or sedimented. Recorded calibration points are indicated by a star (e.g. *point 1) and remain saved even if the device is switched off in the meantime.

1. Fill a black container with a representative sample of the measuring medium. Ensure that it is mixed evenly.
2. Use a part of the sample for the solids analysis in the laboratory.
3. Clean the probe.
4. Dip the probe into the container.
5. Select **MENU**, confirm with **ENTER/ON**.
6. Select **CALIBRATE**, confirm with **ENTER/ON**.
The instrument reads the probe data and info text appears on the display.
7. Select the curve shown and confirm with **ENTER/ON**.
8. Select the desired curve and confirm with **ENTER/ON**.
9. Select **MEMORY** and confirm with **ENTER/ON**.
10. Select **POINT...** (Point 1, 2 or 3), confirm with **ENTER/ON**. In doing so, stir the measuring medium with the probe.

The distance between the probe head and the walls and base of the container must always be more than 70 mm (2.76 in.) (Figure 9).

It takes 5 to 20 seconds to record the calibration point; then the selection menu for calibration appears. While recording, "Memory" is indicated on the display.

11. Take the probe out of the container and clean it. If necessary, record more calibration points with diluted or settled samples. Remove part of these samples for the solids analysis in the laboratory.
12. Determine the solids content of the sample(s) in the laboratory.
13. Select ***POINT...** (*Point 1, 2 or 3) and confirm with **ENTER/ON**.

14. Enter the laboratory values using the navigation keys and confirm with **ENTER/ON**.

- Navigation key **LEFT/RIGHT**: Jump to the next/previous decimal place
- Navigation key **UP/DOWN**: Change number

The laboratory value is saved. The main menu is shown.

Repeat steps 1 to 14 to record additional calibration points.

The instrument automatically sorts the saved calibration points according to the size of the calibration values. The sequence in which the calibration points were recorded is not considered.

Point 1 is always allocated the lowest calibration value.

Point 2 is allocated the next larger calibration value.

Point 3 is allocated the largest calibration value.

The value determined in the laboratory can be corrected at any time by overwriting (refer to [section 5.1.3](#)).

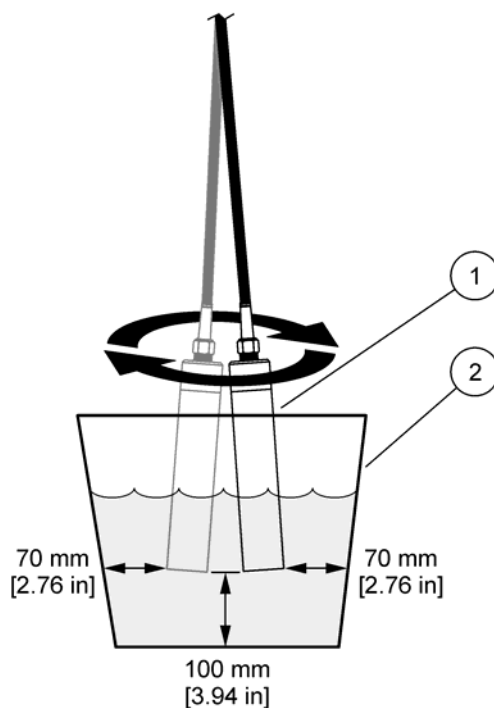


Figure 9 Minimum distances of the probe head in the container

1	Probe	2	Container
---	-------	---	-----------

5.1.3 Manual calibration value correction

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **CALIBRATE**, confirm with **ENTER/ON**.
The instrument reads the probe data.
3. Select the curve shown and confirm with **ENTER/ON**.
4. Select the desired curve and confirm with **ENTER/ON**.
5. Select the desired point and confirm with **ENTER/ON**.
6. Overwrite the existing calibration point.
7. Enter the value for the calibration point using the navigation keys and confirm with **ENTER/ON**.
 - Navigation key **LEFT/RIGHT**: Jump to next/previous decimal place
 - Navigation key **UP/DOWN**: Change number

5.1.4 Delete calibration point

Same procedure as in [section 5.1.3](#). Set the concentration of the calibration point to the value 0.000 and confirm with **ENTER/ON**. The calibration point is deleted.

5.1.5 Reset to default calibration

If the turbidity calibration was changed, it can be reset to the default calibration in accordance with ISO 7027.

Note: Stop the current measurement before the default calibration is reset.

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **CALIBRATE**, confirm with **ENTER/ON**.
The instrument reads the probe data.
3. Select turbidity curve **C-Tu** and confirm with **ENTER/ON**.
Note: If turbidity curve C-Tu is not shown in the menu, select the curve shown and confirm with ENTER/ON. The selection menu for the curves is shown. Select curve C-Tu and confirm with ENTER/ON. Curve C-TU is activated.
Open the selection menu for curve C-TU with ENTER/ON, select RESET and confirm with ENTER/ON.

4. Press navigation key **UP/DOWN**, select **RESET** and confirm with **ENTER/ON**.

The DS curves can also be reset to the default calibration (1-point calibration at 4 g/L). However, a precise determination of the solids concentration can only be obtained through calibration (refer to [section 5.1, page 23](#)).

5.2 Measure

Note: In order to obtain precise solids measurement results, you must calibrate the instrument before the first measurement is taken.

Without any calibration of the DS curves, the instrument returns to 1-point default calibration at 4 g/L. Deviations from the actual measurement value may occur.

5.2.1 Selecting the calibration curve

Before starting the measurement, the calibration curve that corresponds to the measurement point must be selected.

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **CALIBRATE**, confirm with **ENTER/ON**.
3. Select the curve shown and confirm with **ENTER/ON**.
4. Select the desired curve with the navigation key **UP/DOWN** and confirm with **ENTER/ON**.
5. Select **CLEAR/OFF** twice to return to the main menu.

5.2.2 Start the measurement

When starting the measurement, the measurement is carried out continuously and the current measurement value is saved at intervals.

1. Select **MEASURE**, confirm with **ENTER/ON**.
2. Select **START** and confirm with **ENTER/ON**.

5.2.3 Interrupt and resume the measurement

1. Press **CLEAR/OFF** for 2 seconds.

The instrument is switched off and the measurement is interrupted.

2. Press **ENTER/ON** for 2 seconds.

The instrument switches on and the measurement continues with the same settings as before the instrument was switched off.

5.2.4 Stop measurement

1. Select **MEASURE**, confirm with **ENTER/ON**.
2. Select **STOP**, confirm with **ENTER/ON**.

5.3 Show the data in the display

The display shows the stored data for the corresponding calibration curves.

1. Select the desired calibration curve (refer to [section 5.2.1, page 27](#)).
2. Select **DATA**, confirm with **ENTER/ON**.
3. Select **DISPLAY DATA** and confirm with **ENTER/ON**.

Only the stored data for the selected calibration curve is displayed.

5.4 Delete data for a calibration curve

1. Select the desired calibration curve (refer to [section 5.2.1, page 27](#)).
2. Select **DATA**, confirm with **ENTER/ON**.
3. Select **DELETE MEP** and confirm with **ENTER/ON**.

Only the stored data for the selected calibration curve is deleted.

4. Select **YES**, confirm with **ENTER/ON**.

Note: **NO** aborts the erasing procedure.

5.5 Delete stored data for all calibration curves

1. Select **MENU**, confirm with **ENTER/ON**.
2. Select **SYSTEM**, confirm with **ENTER/ON**.
3. Enter the value **379** using the navigation keys and confirm with **ENTER/ON**.

⚠ CAUTION

Potential chemical and biological eye and skin hazards.

Only qualified personnel should perform the tasks described in this section of the operating manual.

The cleanliness of the measurement windows in the sensor head is critical to the accuracy of the measurement results!

Cleaning the instrument

Clean the instrument with a damp, lint-free cloth.

Cleaning the measurement windows

The windows are made of sapphire glass. If required, they can be cleaned with any conventional cleaning agent and a soft cloth. In the case of very stubborn deposits, clean with a soft cloth soaked in 5% hydrochloric acid.

Observe the safety conditions and wear protective clothing!

- Safety glasses
- Gloves
- Overalls

Section 7 Troubleshooting

7.1 Error messages

	Caused by	Action
Communication with probe	The probe is not connected to the instrument	<ul style="list-style-type: none"> – Plug in the probe – Check the connection to the probe – Correct probe number
EEPROM error	Values cannot be saved	Switch the instrument off and then back on again
Overflow area	Concentration too high	<ul style="list-style-type: none"> – Probe is measuring in an unsuitable medium – It may be necessary to recalibrate
Clean the probe	Probe is dirty	<ul style="list-style-type: none"> – Clean the probe – May be triggered by inhomogeneous measuring medium
Offset error	Calibration point set incorrectly (calibration point too close to the zero point)	Recalibration required
Calibration error	No calibration available, calibration points set incorrectly (upper and lower calibration points transposed)	Recalibration required. (1-point default calibration of DS curves not sufficiently precise for measuring medium)
Calibration insufficient –	Further calibration point with lower concentration is missing	Further calibration point with lower concentration is required
Calibration insufficient +	Further calibration point with greater concentration is missing	Further calibration point with greater concentration required
Probe error	Technical error in the probe	<ul style="list-style-type: none"> – Switch the instrument off and then back on again – If the message is repeated, send the probe for repair
Battery is empty	Battery not charged	Charge the battery fully
Incorrect probe	Probe has been replaced but not re-registered or correctly registered on instrument	<ul style="list-style-type: none"> – Connect the correct probe – Use code 390 to release DS curve for new probe and recalibrate
Error time/date	Time/date implausible	<ul style="list-style-type: none"> – Switch the instrument off and then back on again – Check the time/date set – If required, change the clock batteries¹⁾

¹ If the message is repeated, send the instrument for repair

7.2 Informative messages

	Caused by	Measure/note
Wait	Calibration is being recorded	Wait (probe in calibration medium)
Manufacturer standard	Measurement value is outside the standard (only for turbidity)	Value is above the DIN standard of 4000 FNU (NTU)

Section 8 Replacement parts and accessories

8.1 Replacement Parts

Description	Order-No.
TSS Portable handheld instrument (TSS Portable handheld measurement instrument includes meter, TSS probe (10 m (32.8 ft) cable, plug), charger with four adapters (ES, USA, GB and AUS/China), rechargeable batteries, manual and carry case	LXV322.99.00001
TSS Portable handheld instrument	LXV320.99.00001
TSS probe (10 m (32.8 ft) cable, plug)	LXV321.99.00001
6x rechargeable NiMH batteries, AA min. 1.8 Ah	LZY604
Charger for rechargeable batteries with four adapters	LZY607
Battery holder	LZY606
Standard case	LZY605

Section 9 Contact information

HACH Company World Headquarters

P.O. Box 389
Loveland, Colorado
80539-0389 U.S.A.
Tel (800) 227-HACH
(800) -227-4224
(U.S.A. only)
Fax (970) 669-2932
orders@hach.com
www.hach.com

HACH LANGE GMBH

Willstätterstraße 11
D-40549 Düsseldorf
Tel. +49 (0)2 11 52 88-320
Fax +49 (0)2 11 52 88-210
info@hach-lange.de
www.hach-lange.de

HACH LANGE GMBH

Rorschacherstrasse 30a
CH-9424 Rheineck
Tel. +41 (0)848 55 66 99
Fax +41 (0)71 886 91 66
info@hach-lange.ch
www.hach-lange.ch

HACH LANGE APS

Åkandevej 21
DK-2700 Brønshøj
Tel. +45 36 77 29 11
Fax +45 36 77 49 11
info@hach-lange.dk
www.hach-lange.dk

HACH LANGE LDA

Av. do Forte nº8
Fracção M
P-2790-072 Carnaxide
Tel. +351 214 253 420
Fax +351 214 253 429
info@hach-lange.pt
www.hach-lange.pt

HACH LANGE KFT.

Vöröskereszt utca. 8-10.
H-1222 Budapest XXII. ker.
Tel. +36 1 225 7783
Fax +36 1 225 7784
info@hach-lange.hu
www.hach-lange.hu

Repair Service in the United States:

HACH Company
Ames Service
100 Dayton Avenue
Ames, Iowa 50010
Tel (800) 227-4224
(U.S.A. only)
Fax (515) 232-3835

HACH LANGE LTD

Pacific Way
Salford
GB-Manchester, M50 1DL
Tel. +44 (0)161 872 14 87
Fax +44 (0)161 848 73 24
info@hach-lange.co.uk
www.hach-lange.co.uk

HACH LANGE FRANCE S.A.S.

8, mail Barthélémy Thimonnier
Lognes
F-77437 Marne-La-Vallée
cedex 2
Tél. +33 (0) 820 20 14 14
Fax +33 (0)1 69 67 34 99
info@hach-lange.fr
www.hach-lange.fr

HACH LANGE AB

Vinthusdsvägen 159A
SE-128 62 Sköndal
Tel. +46 (0)8 7 98 05 00
Fax +46 (0)8 7 98 05 30
info@hach-lange.se
www.hach-lange.se

HACH LANGE SP. ZO.O.

ul. Krakowska 119
PL-50-428 Wrocław
Tel. +48 801 022 442
Zamówienia: +48 717 177 707
Doradztwo: +48 717 177 777
Fax +48 717 177 778
info@hach-lange.pl
www.hach-lange.pl

HACH LANGE S.R.L.

Str. Căminului nr. 3,
et. 1, ap. 1, Sector 2
RO-021741 București
Tel. +40 (0) 21 205 30 03
Fax +40 (0) 21 205 30 17
info@hach-lange.ro
www.hach-lange.ro

Repair Service in Canada:

Hach Sales & Service
Canada Ltd.
1313 Border Street, Unit 34
Winnipeg, Manitoba
R3H 0X4
Tel (800) 665-7635
(Canada only)
Tel (204) 632-5598
Fax (204) 694-5134
canada@hach.com

HACH LANGE LTD

Unit 1, Chestnut Road
Western Industrial Estate
IRL-Dublin 12
Tel. +353(0)1 460 2522
Fax +353(0)1 450 9337
info@hach-lange.ie
www.hach-lange.ie

HACH LANGE NV/SA

Motstraat 54
B-2800 Mechelen
Tel. +32 (0)15 42 35 00
Fax +32 (0)15 41 61 20
info@hach-lange.be
www.hach-lange.be

HACH LANGE S.R.L.

Via Rossini, 1/A
I-20020 Lainate (MI)
Tel. +39 02 93 575 400
Fax +39 02 93 575 401
info@hach-lange.it
www.hach-lange.it

HACH LANGE S.R.O.

Zastrčená 1278/8
CZ-141 00 Praha 4 - Chodov
Tel. +420 272 12 45 45
Fax +420 272 12 45 46
info@hach-lange.cz
www.hach-lange.cz

HACH LANGE

8, Kr. Sarafov str.
BG-1164 Sofia
Tel. +359 (0)2 963 44 54
Fax +359 (0)2 866 15 26
info@hach-lange.bg
www.hach-lange.bg

Repair Service in Latin America, the Caribbean, the Far East, Indian Subcontinent, Africa, Europe, or the Middle East:

Hach Company World
Headquarters,
P.O. Box 389
Loveland, Colorado,
80539-0389 U.S.A.
Tel +001 (970) 669-3050
Fax +001 (970) 669-2932
intl@hach.com

HACH LANGE GMBH

Hütteldorfer Str. 299/Top 6
A-1140 Wien
Tel. +43 (0)1 912 16 92
Fax +43 (0)1 912 16 92-99
info@hach-lange.at
www.hach-lange.at

DR. LANGE NEDERLAND B.V.

Laan van Westroijen 2a
NL-4003 AZ Tiel
Tel. +31(0)344 63 11 30
Fax +31(0)344 63 11 50
info@hach-lange.nl
www.hach-lange.nl

HACH LANGE S.L.U.

Edificio Seminario
C/Larrauri, 1C- 2ª Pl.
E-48160 Derio/Vizcaya
Tel. +34 94 657 33 88
Fax +34 94 657 33 97
info@hach-lange.es
www.hach-lange.es

HACH LANGE S.R.O.

Roľnícka 21
SK-831 07 Bratislava –
Vajnory
Tel. +421 (0)2 4820 9091
Fax +421 (0)2 4820 9093
info@hach-lange.sk
www.hach-lange.sk

HACH LANGE SU ANALIZ SİSTEMLERİ LTD.ŞTİ.

İlkbahar mah. Galip Erdem
Cad. 616 Sok. No:9
TR-Oran-Çankaya/ANKARA
Tel. +90312 490 83 00
Fax +90312 491 99 03
bilgi@hach-lange.com.tr
www.hach-lange.com.tr

Contact information

HACH LANGE D.O.O.

Fajfarjeva 15
SI-1230 Domžale
Tel. +386 (0)59 051 000
Fax +386 (0)59 051 010
info@hach-lange.si
www.hach-lange.si

HACH LANGE E.Π.E.

Αυλίδος 27
GR-115 27 Αθήνα
Τηλ. +30 210 7777038
Fax +30 210 7777976
info@hach-lange.gr
www.hach-lange.gr

HACH LANGE D.O.O.

Ivana Severa bb
HR-42 000 Varaždin
Tel. +385 (0) 42 305 086
Fax +385 (0) 42 305 087
info@hach-lange.hr
www.hach-lange.hr

HACH LANGE MAROC SARLAU

Villa 14 – Rue 2 Casa
Plaisance
Quartier Racine Extension
MA-Casablanca 20000
Tél. +212 (0)522 97 95 75
Fax +212 (0)522 36 89 34
info-maroc@hach-lange.com
www.hach-lange.ma

HACH LANGE OOO

Finlyandsky prospekt, 4A
Business Zentrum "Petrovsky
fort", R.803
RU-194044, Sankt-Petersburg
Tel. +7 (812) 458 56 00
Fax. +7 (812) 458 56 00
info.russia@hach-lange.com
www.hach-lange.com

Section 10 Limited warranty

Hach Company warrants its products to the original purchaser against any defects that are due to faulty material or workmanship for a period of one year from date of shipment unless otherwise noted in the product manual.

In the event that a defect is discovered during the warranty period, Hach Company agrees that, at its discretion, it will repair or replace the defective product or refund the purchase price excluding original shipping and handling charges. Any product repaired or replaced under this warranty will be warranted only for the remainder of the original product warranty period.

This warranty does not apply to consumable products such as chemical reagents; or consumable components of a product, such as, but not limited to, lamps and tubing.

Contact Hach Company or your distributor to initiate warranty support. Products may not be returned without authorization from the Hach Company.

Limitations

This warranty does not cover:

- Damage caused by acts of God, natural disasters, labor unrest, acts of war (declared or undeclared), terrorism, civil strife or acts of any governmental jurisdiction
- Damage caused by misuse, neglect, accident or improper application or installation
- Damage caused by any repair or attempted repair not authorized by the Hach Company
- Any product not used in accordance with the instructions furnished by the Hach Company
- Freight charges to return merchandise to the Hach Company
- Freight charges on expedited or express shipment of warranted parts or products
- Travel fees associated with on-site warranty repair

This warranty contains the sole express warranty made by the Hach Company in connection with its products. All implied warranties, including without limitation, the warranties of merchantability and fitness for a particular purpose, are expressly disclaimed.

Some states within the United States do not allow the disclaimer of implied warranties and if this is true in your state the above limitation may not apply to you. This warranty gives you specific rights, and you may also have other rights that vary from state to state.

This warranty constitutes the final, complete, and exclusive statement of warranty terms and no person is authorized to make any other warranties or representations on behalf of Hach Company.

Limitation of Remedies

The remedies of repair, replacement or refund of purchase price as stated above are the exclusive remedies for the breach of this warranty. On the basis of strict liability or under any other legal theory, in no event shall the Hach Company be liable for any incidental or consequential damages of any kind for breach of warranty or negligence.

